


PROTECT[®] AND CARE

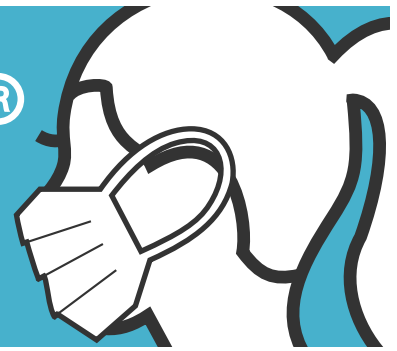

PRODUCTFICHE

Medische nitril handschoenen 3.5g


Non steriel

Ref. A00136

Specificaties

- Voor éénmalig gebruik
- Geschikt voor HACCP en medische doeleinden
- Medisch & beschermend: voorkomen besmettingsgevaar en beschermen tegen virussen, bacteriën & fungi
- 100% Nitril: sterker dan andere handschoen- materialen zoals latex, natuur rubber of vinyl.
- Zijn flexibel en bieden een veel grotere tastgevoeligheid dan standaard diagnostische handschoenen.
- Goede punctie- en scheurbestendigheid
- Hypo allergeen: - Latexvrij
- Poedervrij: minder kans op irritatie.
- Ondoorlatend bij gebruik in natte of vette werkomgeving
- Comfortabele pasvorm:
 - voorzien van rolrand voor makkelijk aan- en uittrekken
 - zowel links als rechts te dragen
- Microgestructureerde vingertoppen

Maten:	Lengte mm	Palmbreedte mm
S 6-7	240	80 +/-10
M 7-8	240	95 +/-10
L 8-9	240	110 +/-10
XL 9-10	240	+110 +/-10

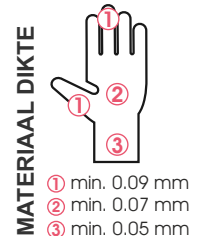

- Duidelijke maataanduiding met kleur
- Verpakking voorzien van barcode en QR code
- Info in 27 talen: NL, FR, EN, DE, ES, IT, PT, BG, CS, DA, EL, FI, HU, HR, KK, LT, LV, NO, PL, RO, RU, SK, SL, SR, BS, SV, TR & UK

Verpakking kwantiteit

Dispenserdoos van 100 stuks/doos 10 doosjes/master box = 1000 handschoenen(500 paar)

EN ISO 374-1:2016
+A1:2018
TYPE B



GESCHIKT VOOR HACCP EN MEDISCHE DOELEINDEN
APPROPRIE POUR HACCP ET SOINS MEDICALES
SUITABLE FOR HACCP AND MEDICAL PURPOSES
FÜR HACCP UND MEDIZINISCHE ZWECKE GEEIGNET



Letter	Chemical	CAS	Permeation Level EN ISO 374-1:2016+A1:2018 Type B	Degradation Level EN ISO 374-4:2019
K	40% Sodium Hydroxide	1310-73-2	6	-8.3%
P	30% Hydrogen Peroxide	7722-84-1	2	34.1%
T	37% Formaldehyde	50-00-0	4	34.3%

Permeation Performance Level	1	2	3	4	5	6
Measured breakthrough time (minutes)	>10	>30	>60	>120	>240	>480

In compliance with: MDR (EU) 2017/745 Class I
PPE Regulation (EU) 2016/425: Category III
EN ISO 21420:2020
EN 455-1-2020, EN 455-2-2015, EN 455-3-2015
ASTM D 6978, ASTM F 1671

EU Type Examination and Module C2 Notified Body 2777
SATRA TECHNOLOGY Europe Limited, Bracetown Business Park,
Cloness Dublin 15, D15 YN2P, Republic of Ireland

Test acc. to ASTM D 6978 - Test Chemotherapy

TEST CHEMOTHERAPY DRUG AND CONCENTRATION	MINIMUM BREAKTHROUGH DETECTION TIME (Specimen 1/23) (Minutes)
Cisplatin (50 mg/ml) (5,000 ppm)	14.7 (15.1, 14.7, 16.8)
Cisplatin 1.0 mg/ml (1,000 ppm)	>240 min.
Cyclophosphamide (Cytosar) 20 mg/ml (20,000 ppm)	>240 min.
Cytarabine 100 mg/ml (100,000 ppm)	>240 min.
Doxorubicin (IDC) 10.0 mg/ml (10,000 ppm)	>240 min.
Doxorubicin Hydrochloride 2.0 mg/ml (2,000 ppm)	>240 min.
Etoposide (Toposar) 20.0 mg/ml (20,000 ppm)	>240 min.
Fluorouracil 50.0 mg/ml (50,000 ppm)	>240 min.
Idarubicin 50.0 mg/ml (50,000 ppm)	>240 min.
Methotrexate 25 mg/ml (25,000 ppm)	>240 min.
Mitomycin C 0.5 mg/ml (500 ppm)	>240 min.
Mitomycin 2.0 mg/ml (2,000 ppm)	>240 min.
Paclitaxel (Taxol) 6.0 mg/ml (6,000 ppm)	>240 min.
Trifluoperazine 10.0 mg/ml (10,000 ppm)	58.8 (110.0, 58.8, 67.0)
Vincristine Sulfate 1.0 mg/ml (1,000 ppm)	>240 min.



ISO 374-5:2016



TESTED IN ACCORDANCE TO:
EN ISO 374-5:2016 and EN ISO 374-4:2019
Degradation results indicate the change in puncture resistance
of the gloves after exposure to the challenge chemical
PROTECTION AGAINST BACTERIA, VIRUSES & FUNGI- PASS

SPECIFICATION OF NITRILE DISPOSABLE GLOVES powder-free						
Size Spec		Small	Medium	Large	X-large	Tolerance
Total Length (mm)		240	240	240	240	min
Width of palm (mm)		85	95	105	115	± 5
Thickness (mm)	Finger	0.09	0.09	0.09	0.09	± 0.02
	Palm	0. 07	0.07	0.07	0.07	± 0.02
	Wrist	0.05	0.05	0.05	0.05	± 0.02
Weight per piece (gram)		3.0	3.5	3.9	4.3	± 0.3
T ensile-strength (≧)min		14Mpa	14Mpa	14Mpa	14MPa	min
Elongation (%)min	Before aging:	500	500	500	500	min
	After aging:	400	400	400	400	min